



Invacare[®] | ***Dacapo***

Flexibility & comfort



Invacare **Dacapo**

Invacare **Dacapo** is a broad range of mattresses designed to match customer's requirements within Invacare's homeware and hospital bed ranges. Selecting from high resilience, environmentally friendly, or high density visco-elastic foams within the 'Standard', 'Comfort' and 'Preventive' ranges, **Dacapo** has been developed to offer optimal cushioning and pressure relief without compromising on hygiene and user friendliness.

Invacare **Dacapo Standard**

Dacapo Standard is a standard mattress comprehensively satisfies basic requirements for good cushioning. High quality, high resilient foam provides a resilient & firm supportive base. Fire retardancy conforming to European Standard EN 597 - 1

Dacapo Standard Cover

- 100% cotton
- Washable at 40 degrees



Features and options



10cm of high resilience foam with a density of 38kg/m^3 and a hardness of 140N.



Zipper


Zippered on three sides

Invacare **Dacapo Basic**

Dacapo Basic is the first incontinence mattress in the series. It fully satisfies basic comfort & hygiene requirements for a sleeping surface. The high density, high resilient, foam provides a firm supportive base, strength & durability. Fire retardancy conforming to European Standard EN 597 - 1 and European Standard EN 597 - 2

Dacapo Basic Cover

- 100% Fluid-proof
- Polyurethane coated polyester
- Breathable
- Two way stretch capacity
- Washable at 95 degrees

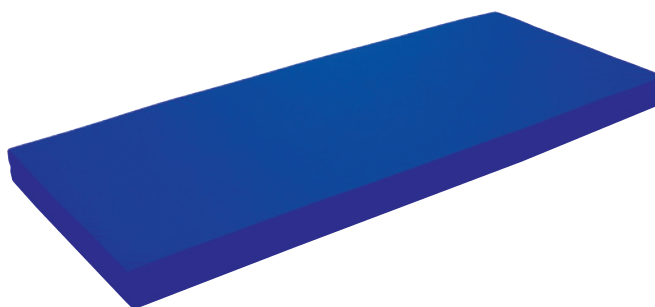

Features and options



12cm of high resilience foam with a density of 38kg/m^3 and a hardness of 140N.



Easy Handling




Zipper

Discreetly zippered on three sides with cloth covering

Invacare **Dacapo Heavy User Light**

Dacapo Heavy User Light is tailored to meet the needs of the larger user and offers a firm, comfortable cushioned support.

Uniquely constructed using two new high resilience foams, this model can be reversed to adjust to the degree of firmness and support required by the user. The broader dimensions enjoyed by **Heavy User Light** make it the perfect complement for Invacare's Broad Home Care Bed models. Fire retardancy conforming to European Standard EN 597 – 1 and European Standard EN 597 – 2

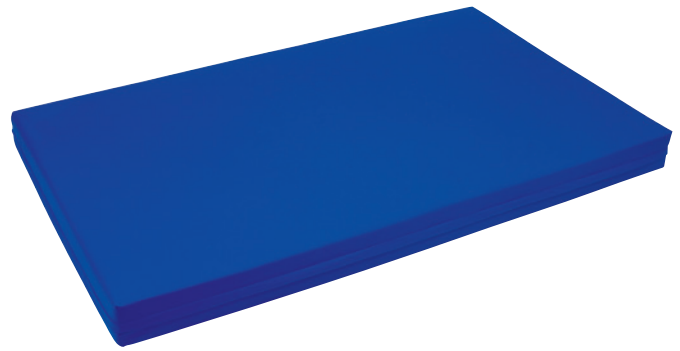


Dacapo Heavy User Light Cover

- 100% fluid-proof
- Polyurethane coated polyester
- Breathable
- Two way stretch capacity
- Washable at 95 degrees

Features and options



Easy Handling

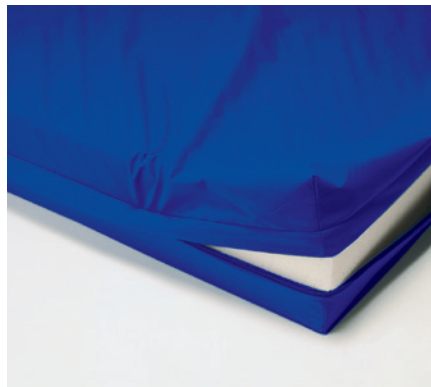

Zipper

Discreetly zippered on three sides with cloth covering



Heavy User Light
maximum 250kg:

Top (5 cm): visco-elastic foam 50kg/m³ density, 100N hardness

Base (5 cm): high resilience foam 38kg/m³ density, 140N hardness

Base (5 cm): high resilience foam 50kg/m³ density, 145N hardness

Heavy User (Available as a special)
maximum 350 kg:

Top (5 cm): visco-elastic foam 50kg/m³ density, 100N hardness

Base (5 cm): high resilience foam 38kg/m³ density, 140N hardness

Base (5 cm): high resilience foam 50kg/m³ density, 200N hardness


Invacare **Dacapo Comfort**

The **Dacapo Comfort** is a comfortable mattress comprised of a foam core whose density provides good recovery while the mattress' unique open cell structure, allows the body to breathe. The mattress is constructed with great flexibility in order to provide a comfortable and supporting surface. The mattress is cut into various sized squares that are individually flexible. This profiling allows the contours of the body to sink into the material so that optimum sleeping comfort is achieved.


Fire retardancy conforming to European Standard EN 597 - 1 and European Standard EN 597 - 2

Dacapo Comfort Cover

- 100% fluid-proof
- Polyurethane coated polyester
- Breathable
- Two way stretch capacity
- Washable at 95 degrees

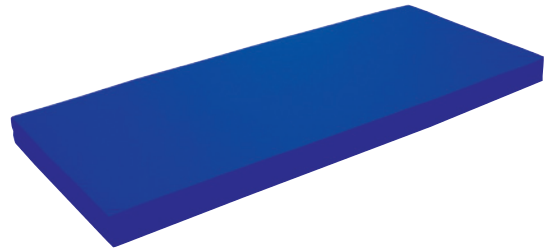

Features and options

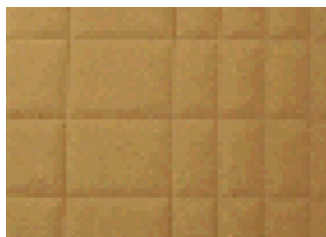

12cm of environmentally friendly foam with a density of 53kg/m³ and a hardness of 145N.



Easy Handling



Zipper

Discreetly zippered on three sides with cloth covering

Invacare **Dacapo Comfort Plus**


The **Dacapo Comfort Plus** is a prevention mattress which satisfies the basic requirements for a comfortable sleep. A unique pressure relieving mattress consisting of various large keyhole cut squares that are individually flexible. The mattress consists of environmentally friendly foam which has a good recovery and its unique open cell structure allows the body to breathe. The mattress offers a dynamic combination of comfort, support and stability while shear and friction forces are efficiently minimized due to the keyhole cut design.

Dacapo Comfort Plus Cover

- 100% fluid-proof
- Polyurethane coated polyester
- Breathable
- Two way stretch capacity
- Washable at 95 degrees



Features and options

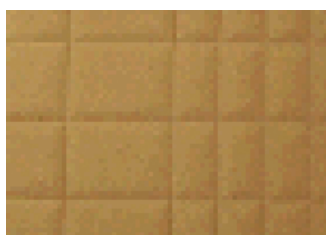

14cm of environmentally friendly foam with a density of 53kg/m³ and a hardness of 145N.



Easy Handling



Zipper

Discreetly zippered on three sides with cloth covering

Invacare **Dacapo Child**

Dacapo Child a single layer mattress constructed from an environmentally-friendly organic foam. The mattress is designed for the smaller, younger user and provides extra cushioning, ultra-durability and enhanced ventilation. In addition, its high volume weight maximises body responsiveness & recovery.


Fire retardency conforming to European Standard EN 597 - 1 and European Standard EN 597 - 2

Dacapo Child Cover

- 100% fluid proof
- Coated polyester with polyurethane*
- Breathable

- Two way stretch capacity
- Washable at 95 degrees

*Also available with cotton cover option.

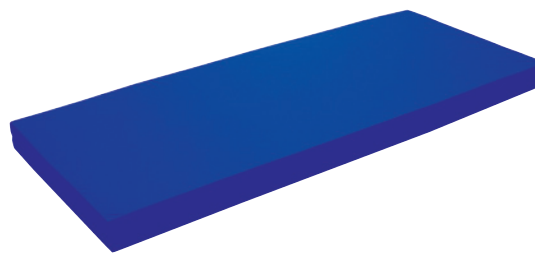

Features and options




Top (8 cm)

Environmentally friendly foam with a density of 50kg/m³ and a hardness of 145N



Easy Handling

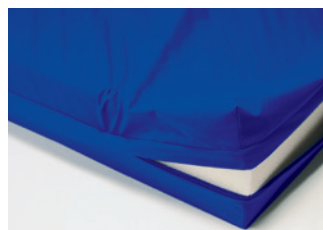


Zipper

Discreetly zippered on three sides with cloth covering

Invacare **Dacapo Square**

Dacapo Square is a combination mattress with an environmentally- friendly, organic foam top. The model's unique square-patterned body contour cut minimises surface pressure, affording the user outstanding pressure relief and body pressure equalisation. The organic foam's high volume weight gives it outstanding responsive & recovery properties, while its unique vapour permeable cell structure guarantees a well-aerated interior. In terms of user benefits, **Dacapo Square** excels in its ability to conform to every curve of the body and provide total body support. Fire retardency conforming to European Standard EN 597- 1 and European Standard EN 597 - 2

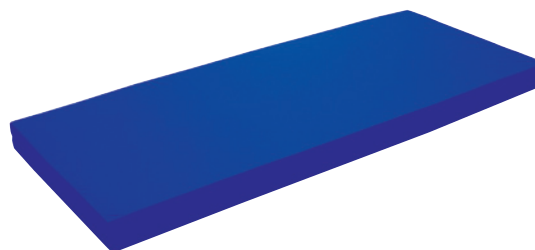


Dacapo Square Cover

- 100% Fluid proof
- Polyurethane coated polyester
- Breathable
- Two way stretch capacity
- Washable at 95 degrees

Features and options




Easy Handling

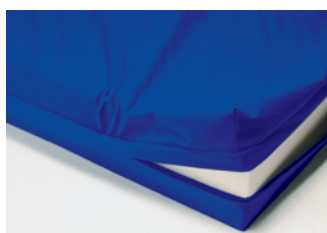

Zipper

Discreetly zippered on three sides with cloth covering

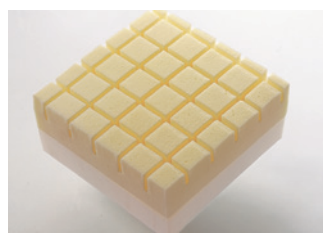

Top (8 cm)

Environmentally friendly foam with a density of 50kg/m³ and a hardness of 145N.

Base (6 cm)

high resilience foam with a density of 38kg/m³ and a hardness of 140N.

Invacare **Dacapo Combi**

Dacapo Combi - is a combination mattress with outstanding pressure relief properties. The top foam is composed of a self-regulating, heat sensitive material which works in conjunction with body temperature & movement. **Dacapo Combi** relaxes according to the user's body dynamics and softens as the user gently sinks into it, comfortably moulding & conforming to the body shape. The open cell structure negates heat & moisture build-up, while the base high resilient, foam reinforces the mattress strength & robustness. Fire retardancy conforming to European Standard EN 597 - 1 and European Standard EN 597 - 2




Dacapo Combi Cover

- 100% fluid proof
- Polyurethane coated polyester
- Breathable
- Two way stretch capacity
- Washable at 95 degrees

Features and options



Easy Handling

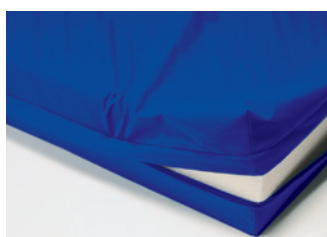

Zipper

Discreetly zippered on three sides with cloth covering



Top (6 cm)


visco-elastic foam with a density of 50kg/m³ and a hardness of 100N.

Base (8 cm)

high resilience foam with a density of 38kg/m³ and a hardness of 140N.

Invacare **Dacapo Top**

Dacapo Top - is a high quality, pressure relieving top layer ideally suited to basic mattresses. **Dacapo Top's** 'temperature smart' foam adapts to the user's body dynamics & temperature. On impact, the material relaxes and then moulds & conforms to the contours of the body leaving the user cradled in comfort. The foam's open cell structure affords a well ventilated interior which reduces the build up of heat & body moisture. Fire retardancy conforming to European Standard EN 597 - 1 and European Standard EN 597 - 2

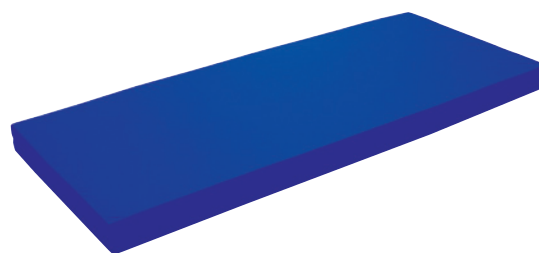

Dacapo Top Cover

- 100% fluid proof
- Polyurethane coated polyester
- Breathable
- Two way stretch
- Washable at 95 degrees

Features and options



Zipper

Discreetly zippered on three sides with cloth covering



6 cm

visco-elastic foam with a density of 50kg/m³ and a hardness of 100N.

Cover Options

Wipedown cover

The wipedown cover is composed of polyutherane coated polyester. The cover is 2 way stretch and gentle on the skin. The fabric is vapour permeable and waterproof. This cover is allergy tested and treated against dust mites.

Evacuation cover

The Evacuation cover is composed of 2 way stretch, vapour permeable, waterproof, polyurethane coated polyester. These covers are equipped with handles on all four sides. On the two long sides are straps, which in case of an emergency evacuation can be wrapped around the patient. This cover is allergy tested and treated against dust mites.


Cotton cover

The Cotton cover is composed of standard cotton and washable at 40 °C.



Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

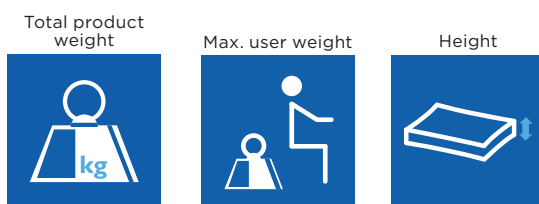

Recommended min and max user weight:

	Total product weight	Max. user weight	Height
Standard	7 - 10.5 kg	40 - 100 kg	100 mm
Basic	8.5 - 12.5 kg	40 - 130 kg	120 mm
Heavy User Light	11 - 16.5 kg	40 - 250 kg	150 mm
Heavy User	12 - 17.5 kg	40 - 350 kg	150 mm
Comfort	11.5 - 17kg	40 - 130 kg	120 mm
Comfort Plus	13.5 - 19kg	40 - 140 kg	140 mm
Child	3 kg	0 - 100 kg	80 mm
Square	10 - 15 kg	40 - 140 kg	140 mm
Combi	8.5 - 13 kg	40 - 140 kg	140 mm
Top	5 - 7.5 kg	40 - 130 kg	60 mm

Frame colours



Scan here to access the website.

Tel: +41 61 487 70 70
hqueurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland



© 2017 Invacare International Sàrl. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Dacapo - EU - 13/2018