



Invacare® | **Medley[®] Ergo**

High quality, cost effective
community profiling bed



The **Medley Ergo** from Invacare offers the perfect balance of economy, functionality and quality. The optimised length of the mattress platform sections provides increased comfort for clients when sitting up in bed and also reduces shear and friction when the bed sections are adjusted. The **Medley Ergo** can be easily transformed into the **Invacare Medley Ergo Select**. This model is particularly stylish with natural beech facias covering the motors, creating a homely feel.



Optimised dimensions of the mattress platform

The mattress platform sections has optimal dimensions designed in accordance with anthropometric statistics and has been pressure mapped during the design process.. The **Medley Ergo** offers:


- Increased comfort when sitting
- Reduced shear and friction from sliding downwards during profiling
- Suitability for shorter clients

Improved functionality and efficiency

The **Medley Ergo** has been specially designed for speedy assembly and quick transportation to meet the needs of Residential and Community care. The **Medley Ergo** can be supplied on a transporter kit so that it can easily be moved and installed. In addition, any assembling can be done without the use of tools.

Height range supports clients and carers

The flexible design of the **Medley Ergo** allows the bed to be mounted at a low height, or a higher position if preferred, creating a range of working heights for the carer and transfer heights for the client. This is facilitated by the unique feature of the bed end, which enables the mattress support to be mounted in both upper and lower positions.

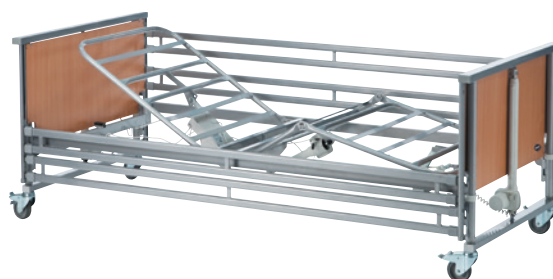

Features and options



Full length side rails

The **Medley Ergo** is available with wooden or steel full length side rails. Wooden extendable side rails (150mm) are also available as an option.



Scala 2 range

A range of collapsible, three quarter length, steel side rails are available in various heights to cover standard to deep mattress heights, such as dynamic mattress systems.

Options



Homely design

The **Medley Ergo Select** offers an attractive, homely design. Available as a retro-fit to the **Medley Ergo**, this aesthetic bed end cover fits with precision over the bed actuators, and can be removed at any time.

The **Medley Ergo Select** has a leg support that is equipped with a rastofix to enable individual positioning for the legs.

Features and Options



Easy to transport

Medley Ergo beds are easy to assemble and dismantle and can be transported using transporter kits (optional).



Lifting Pole

Can be adjusted in both height and depth (optional).



Select design for Medley Ergo

The homely design of **Medley Ergo Select** can also be retro-fitted (optional).



Double height positions for bed ends

Adjustable height range for **Medley Ergo** of 330-730mm or 400-800mm.



Optimised dimensions of mattress platform

According to anthropometrical measurements, the optimised dimensions meets the body size of a wider number of clients and reduces shear and friction.

Features and options



Mattress support extension

Extension kit for increasing the length of the bed by 150mm.



Transport fitting

For easy storage and transport.















Hand controls

Available with or without
lockable functions.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk.

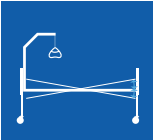Width outside	Width inside	Length outside	Length inside	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle	Legrest lift
								
1020 mm	900 mm	2240 mm	2080 mm	400 - 800 mm or 330 - 730 mm	800 - 230 - 300 - 670 mm	0 - 70°	0 - 24°	0 - 11°

				Complies to the safety standard, IEC 60601-2-52.
Tilt option ¹	Total weight	Max. user weight	Heaviest part of the product	
11°	64 kg	SWL 180 kg 145 kg	17 kg	

1. Please note: the tilt function is only used for clients under medical supervision. If there is a risk of clients using the tilt function without supervision, please use a lockable hand control.

NB

The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

Light grey

Beech



Please note colours may vary slightly from those displayed above.

Warranty

3 years frame
3 years actuator
2 years handset

The *Invacare* **Medley Ergo** is part of the Safe Patient Handling range.



Safe Patient Handling

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

